


# TPWsoft

SOTAX

Automated Sample Preparation allows for robust and reproducible preparations that are equivalent to or better than your manual procedures. Secure sample tracking in compliant software tracks and records every step in your automated workflow.

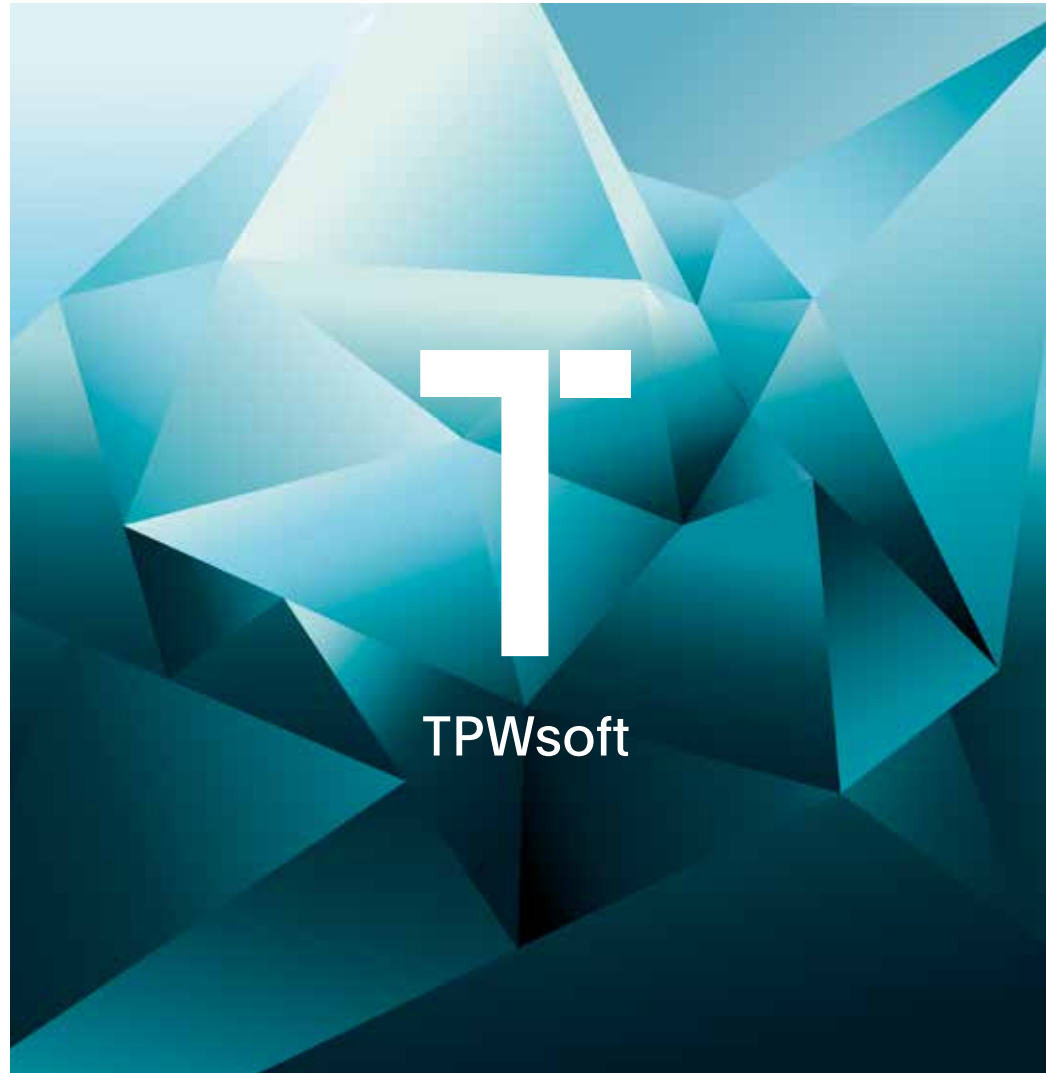

## Sample Preparation Software



# Automate your method development.

The user-friendly software interface allows for rapid adoption in all environments from academic AR&D to the QC lab. Advanced developer options accelerate the method development process to facilitate Analytical Quality by Design (AQbD).

The software architecture with central MS-SQL database is 100% scalable and fulfills all requirements for implementation of a 21 CFR part 11 compliant system. The intuitive software interface streamlines the method transfer across sites. EmpowerLink™ can be used to automatically link sample information to valid results with native Empower sample set integration in both online and offline HPLC workflows.

100% scalable software architecture

21 CFR part 11 compliant

Easy-to-use & intuitive software interface

EmpowerLink™



sotax.com/apwsoft



sotax.com/tpwsoft

T




TPWsoft

A

APWsoft

T


TPWsoft

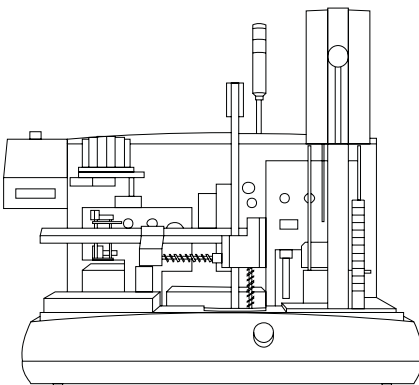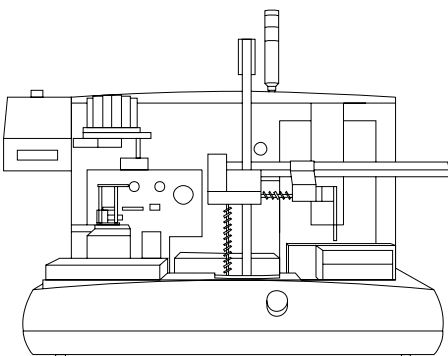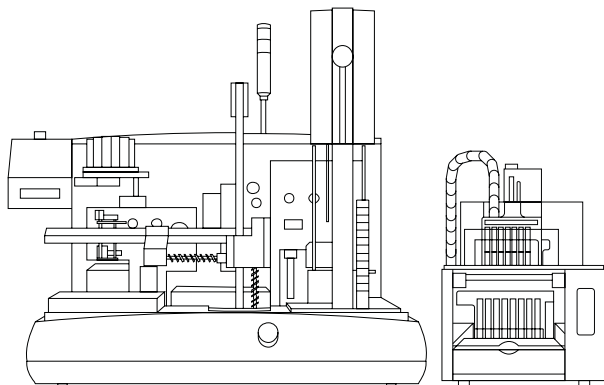




# Framework. Modules. Drivers.

Example: APWsoft & TPWsoft with 3 automated systems and a shared database

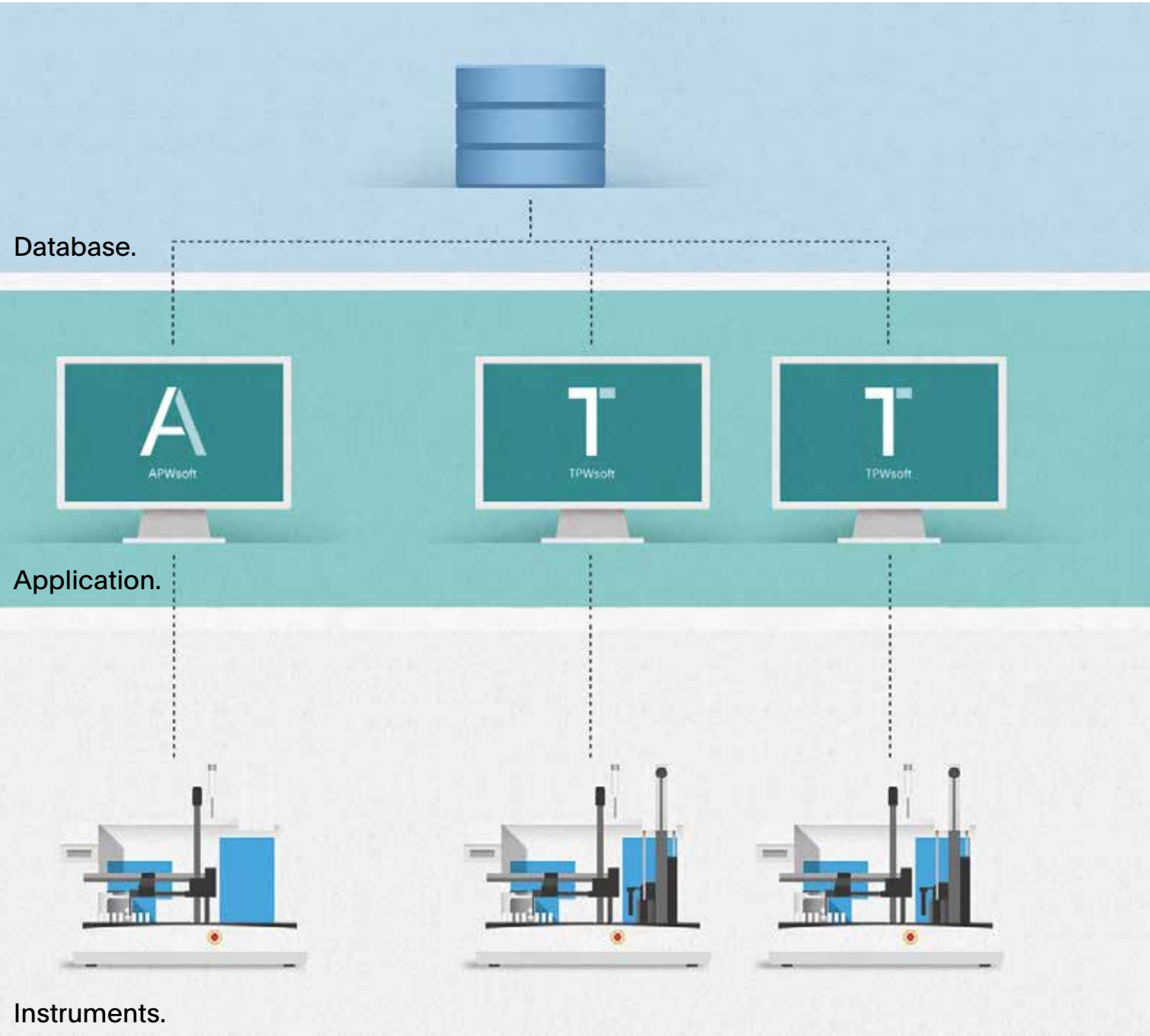

## Framework — All key functions.



Whether installed on a single PC or on multiple networked computers, the APWsoft / TPWsoft framework with central SQL database is designed for both stand-alone and networked operation. It includes all basic functions to collect, record, analyze, evaluate, report, and manage test data in full compliance with 21 CFR part 11 requirements.

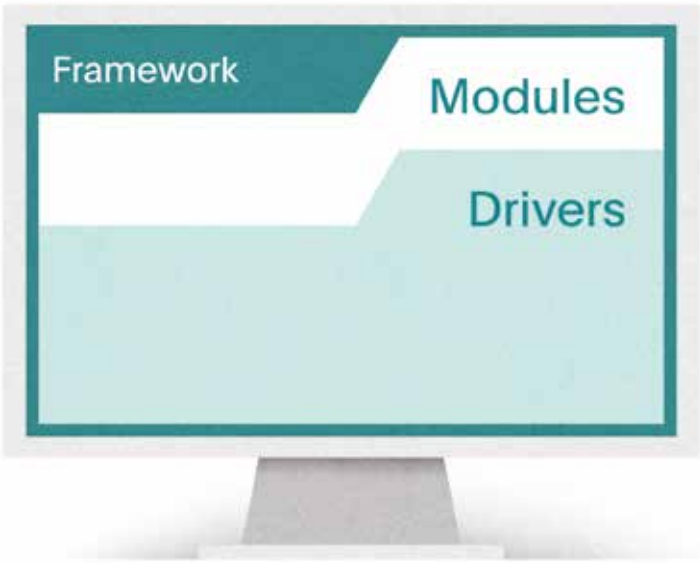



# EmpowerLink™ — Online.

The automated system can be directly connected to an HPLC or Waters Acquity UPLC. After each sample preparation, the online system will prepare and inject fresh samples. When using the Waters Empower Chromatography Data System (CDS), the sample information can be sent directly to the CDS to automatically create a sample set during sample injection. Important information during preparation (such as ID, weight, and sample concentration) can be automatically sent to CDS for analysis.

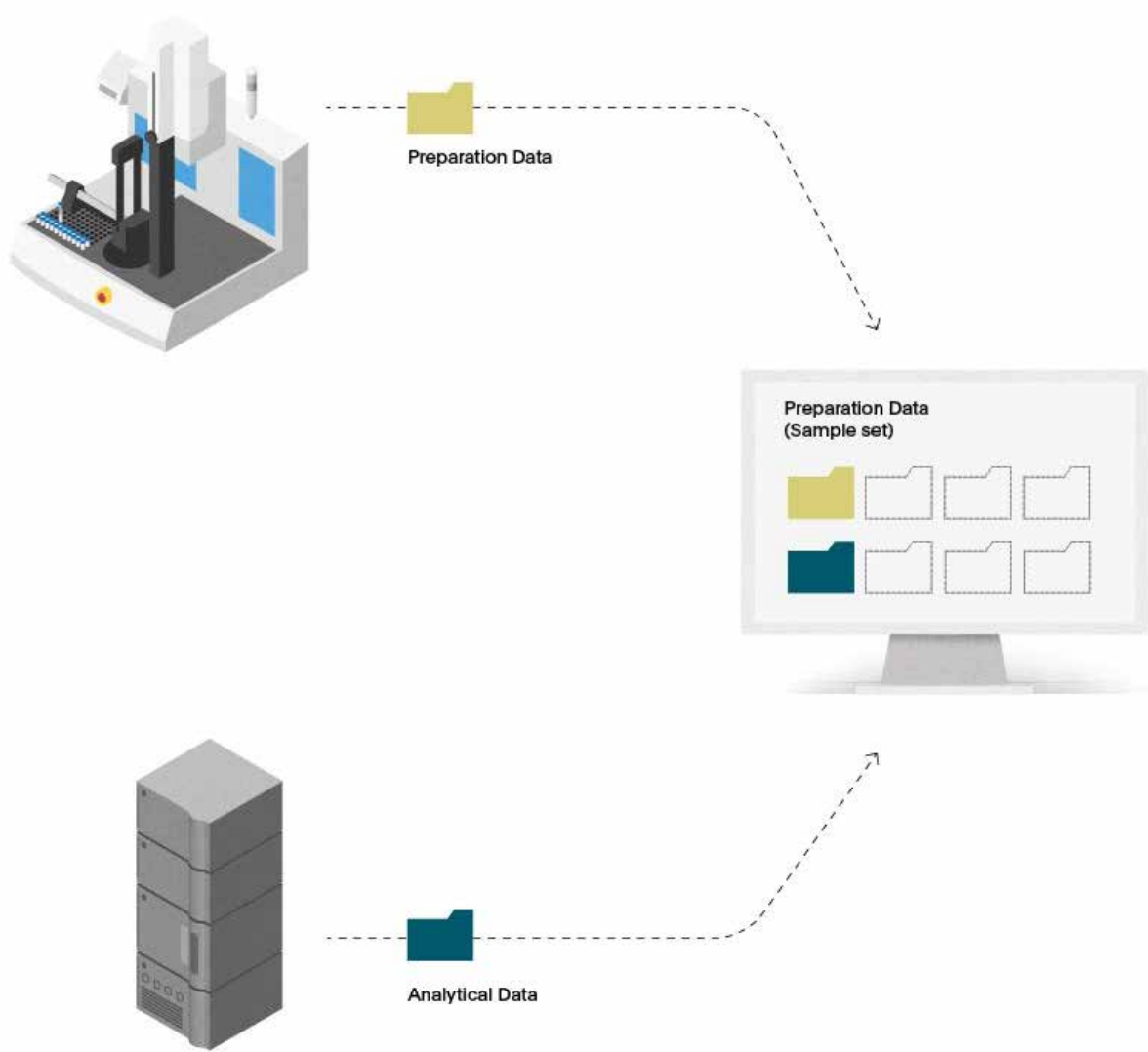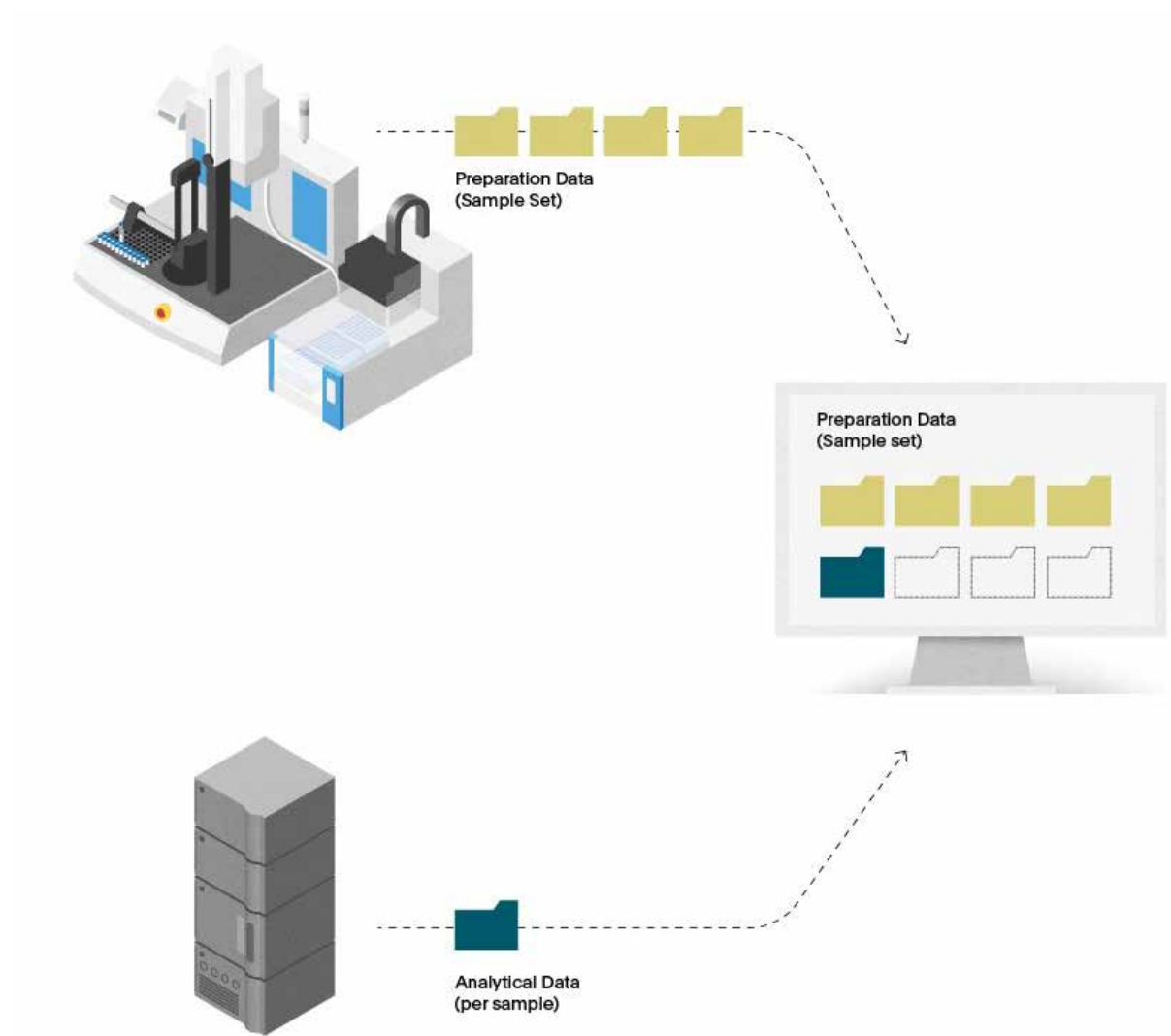

# EmpowerLink™ — Offline.

Empower sample sets can be created automatically without a physical instrument connection. The samples can be collected offline in the LC vial plates on the SAM Sample Manager and later moved to the LC autosampler. EmpowerLink Offline can be used to create sample sets on the HPLC and transmit all relevant sample information with zero vial re-racking.

APWsoft & TPWsoft – Technical Specifications

TPWsoft / APWsoft		
Supported Windows versions		
Windows 10 64-bit only		
Controls (minimum requirements)	CPU	Intel Core i5, 8th Gen or better
	RAM	8 GB or more
	Hard disk space	At least 10 GB
	Video resolution	Minimum resolution of 1920 × 1080 pixel (HD resolution)
	Ethernet	At least 100 Mbit Ethernet connection or better
Supported printers		
All printers supported by OS		
Database		
MS SLQ Server 2014 or greater		
Accessories: Cameras		
2× independent channel USB 2.0 slots		

Technical specifications are subject to change without prior notice. Products illustrated in this brochure may include options or modifications not fitted as standard. No liability for errors and omissions.

SOTAX Worldwide.

Europe

Switzerland (HQ Europe)  
Aesch / Basel  
P Tech Support +41 61 487 5460  
P Office +41 61 487 5454  
info@sotax.com

Czech Republic  
Prague  
P Tech Support +41 61 487 5460  
P Office +420 246 039 260  
sotaxcz@sotax.com

France  
Saint-Louis  
P Tech Support +41 61 487 5460  
P Office +33 3 8970 0846  
info@sotax.com

Germany  
Lörrach  
P Tech Support +41 61 487 5460  
P Office +49 7621 16 5635  
info@sotax.com

Great Britain  
Foston / Derbyshire  
P Tech Support +44 1332 981 571  
P Office +44 1332 981 303  
sotaxgb@sotax.com

Italy  
Milan  
P Tech Support +41 61 487 5460  
P Office +39 02 8363 2641  
sotaxitalia@sotax.com

Ireland  
Longford  
P Tech Support +353 43 334 7779  
P Office +353 43 334 7779  
sotaxie@sotax.com

Americas

USA (HQ Americas)  
Westborough, MA  
P Tech Support +1 508 544 4040  
P Office +1 508 417 1112  
sotaxusa@sotax.com

Canada  
Westborough, MA  
P Tech Support +1 508 544 4040  
P Office +1 508 417 1112  
sotaxcanada@sotax.com

Asia-Pacific

Thailand (HQ Asia-Pacific)  
Bangkok  
P Tech Support +91 99 202 11 211  
P Office +66 20 02 71 15  
sotaxasia@sotax.com

India  
Mumbai  
P Tech Support +91 99 202 11 211  
P Office +91 22 4295 0191  
sotaxindia@sotax.com

